

## How to Use Radara

1. Cleanse skin thoroughly.
2. Remove a single pack of Radara™ patches from box and tear off the logo panel. Notice how the correct-shaped patch will be used for each side of the face as presented, using one side at a time.
3. Place a patch on the outer eye area. Once in position, **firmly smooth the patch** over the skin with finger pressure for five seconds - you should feel a painless sensation.
4. Replace each patch into the pack, needle side up.
5. Dispense a single pump of the serum onto a clean finger and firmly smooth over the prepared skin, focusing on the crow's feet/patch area and being careful to not get any in the eye.
6. Repeat step 3 and now leave in place for five minutes. The patches should now hold in place on the skin, however try to minimise facial movements.
7. After five minutes, remove and dispose of the patches and wipe away any excess serum. **Do not re-use patches.**

The Radara™ skincare regime consists of daily patches applied for a 4-week period and is recommended for use three to four times per year for maximum results.

*NOTE: Avoid contact with the eyes. If product gets in the eyes, rinse well with water immediately. It is common to observe a mild redness when using the patches - this will subside in less than 15 minutes. On first use the pump has to be primed once only.*

Radara™ is manufactured in the UK by Innoture Medical Technology Ltd. For more information and video instructions see [www.radara.co.uk](http://www.radara.co.uk)



Radara™ has been designed, developed and manufactured by Innoture Medical Technology Ltd, a company based in Swansea, Wales, UK, with research directed by leading skincare advisors based at the state-of-the-art Institute of Life Science at Swansea University.

Contact Us  
Innoture Medical Technology  
15a Ives St, London, SW3 2ND  
+44 (0)207 594 4947  
[info@radara.co.uk](mailto:info@radara.co.uk)

Visit us on the web  
[www.radara.co.uk](http://www.radara.co.uk)  
[@RadaraUK](https://www.instagram.com/RadaraUK)



^  
RADARA  
Targeted Skincare

Rejuvenating Micro-Channelling  
Patches to Gently Correct  
Fine Lines & Wrinkles

PATIENT INFORMATION LEAFLET

## Radara


*noun*

Combination concept of 'R' for radiance and 'Adara' - from the Greek meaning 'beautiful'

### Targeted, Sophisticated Skincare


- The Radara™ skincare system consists of rejuvenating micro-channelling patches and high purity hyaluronic acid (HA) serum to restore skin radiance and quality, whilst diminishing signs of ageing
- This innovative, targeted and painless approach delivers smoother, firmer skin with reduced lines and wrinkles in just four weeks
- Radara™ is designed for once-daily use from the comfort of your own home
- The treatment regime lasts for four weeks, with the benefits being seen for up to eight weeks in total
- Radara™ works by creating micro-channels through the upper layers of the skin, to bypass the skin barrier and allow HA to penetrate to the deeper levels where it is most needed to restore skin's natural elasticity, hydration and support

### How Does Radara Work?

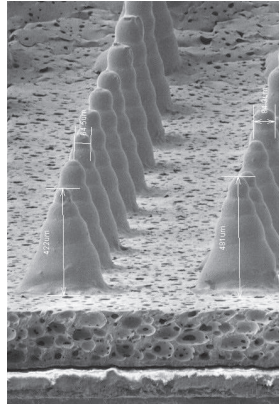

The cleverly shaped, flexible Radara patches are coated with microscopic plastic structures (like needles) less than 0.5mm long. When applied to the skin, these patches painlessly create hundreds of tiny micro-channels in the epidermis.

Radara's specially-formulated serum containing pure hyaluronic acid (HA) is then applied and flows through the newly-created micro-channels to the deeper skin layers where it is needed most. Here, the HA serum restores the skin's natural elasticity, hydration and support – giving a smoother, firmer feel as well as gently reducing signs of ageing such as fine lines and wrinkles.

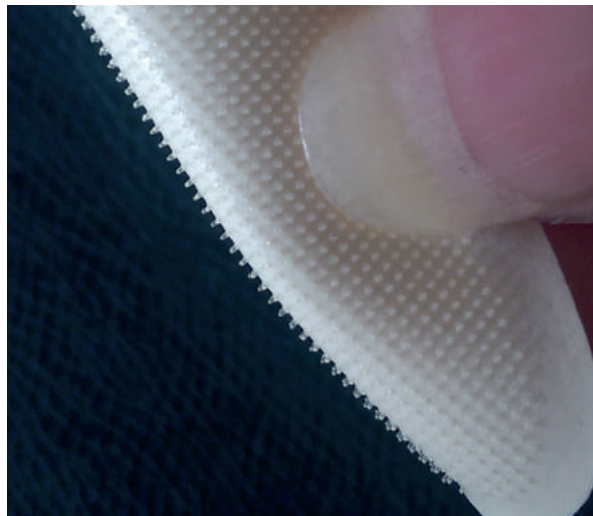

Radara™ works deep within the skin, bringing real changes from the inside out by directly improving skin cell renewal and hydration.

### Results

**An 8-week trial independently assessed by clinical dermatologists showed:**

- Significant reduction in fine lines and wrinkles of up to 35% in four weeks, with noticeable improvements in just two weeks
- Radara™ patches and HA serum were almost twice as effective at reducing fine lines and wrinkles compared to using serum alone
- Radara™ continued to deliver results for up to four weeks post-treatment (improved skin quality, wrinkle reduction)
- Five-minute daily application with no downtime or adverse side effects - unlike traditional microneedling, which may require pre-treatment anaesthetic and can lead to bleeding

### Testimonials

*"I love the results I've seen with Radara!"*

*"I wanted something gentle but effective to give me a lift. Radara was so easy to use!"*

*"There was no pain or discomfort at all, and in just a few weeks, my skin was brighter, fresher and firmer."*

*"I still looked like me, just with more of a youthful glow."*

